

Designed & manufactured to ensure a minimally invasive yet consistently effective procedure.

FEATURES

- Smooth tip reduces risk of trauma upon insertion and removal of catheter.
- Stable placement sheath decreases the need for a tenaculum.
- Soft, concentric balloon secures above the internal os preventing leakage of sterile saline or water-based contrast media.
- Shapeable HS catheter is equipped with an embedded stylet in cases when stenosis is present or when the vaginal tract and cervical canal are not aligned.
- Made with 20% Radiopaque material allowing catheter to be seen under xray.

Download IFU



TECHNICAL DATA

French Size 5Fr, 7Fr
 Length 28cm
 Maximum Balloon Volume 5Fr, 1.5ml/7Fr, 3ml
 Box Qty. 10 ea.
 Sterilization Method R (Gamma)
 Standard Procedures Hysterosalpingography, Saline Infusion Sonography, Sonohysterography

ORDERING INFORMATION

TMI1154 5Fr Flexible HS Catheter
 TMI1155 7Fr Flexible HS Catheter
 TMI1185 5Fr Shapeable HS Catheter
 TMI1187 7Fr Shapeable HS Catheter

All Thomas Medical HS Catheters are PVC Free and Non-DEHP. Not made with natural rubber latex.

SPECIALTY HS CATHETERS



STEEN OPEN-TIP HS CATHETER

The Steen Open-Tip HS catheter utilizes a balloon closer to the distal tip in order to provide a clearer view of the uterus while performing a hysterosalpingography or sonohysterography procedure. The open distal end allows a forward injection into the uterus.

ORDERING INFORMATION

TMI1188 5Fr Steen Open-Tip HS Catheter
 TMI1189 7Fr Steen Open-Tip HS Catheter




MILLER ADVANCE HS CATHETER

The Miller Advance Catheter was designed with features for improved patient comfort with a soft catheter tip and slender catheter shaft that reduces the need for dilation.

ORDERING INFORMATION

TMI2000 5Fr, 21 cm Miller Advance HS Catheter
 TMI2001 5Fr, 25cm Miller Advance HS Catheter
 TMI2002 5Fr, 28cm Miller Advance HS Catheter

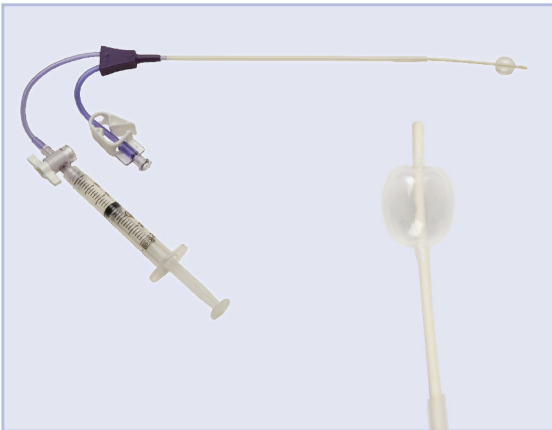


TWO-TAP HS CATHETER

The Two-Tap HS catheter utilizes a second stopcock instead of a pinch clamp for easier manipulation when releasing the saline or water-based contrast media from the uterus.

ORDERING INFORMATION

TMI1181 5Fr Two-Tap HS Catheter
 TMI1182 7Fr Two-Tap HS Catheter
 TMI1183 5Fr Shapeable Two-Tap HS Catheter
 TMI1184 7Fr Shapeable Two-Tap HS Catheter

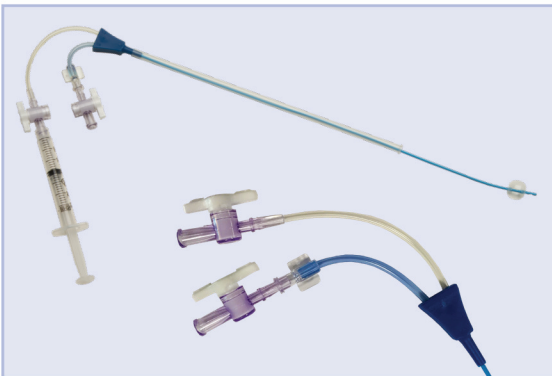

TECHNICAL DATA

CATHETER TYPE	FRENCH SIZE(S)	LENGTH (CM)	MAXIMUM BALLOON VOLUME (ml)	BOX QTY	STERILIZATION METHOD
FLEXIBLE HS CATHETER	5FR/7FR	28	5Fr, 1.5ml/7Fr, 3ml	10	R (GAMMA)
SHAPEABLE HS CATHETER	5FR/7FR	28	5Fr, 1.5ml/7Fr, 3ml	10	R (GAMMA)
STEEN OPEN-TIP HS CATHETER	5FR/7FR	28	5Fr, 1.5ml/7Fr, 3ml	10	R (GAMMA)
MILLER ADVANCE HS CATHETER	5FR	21, 25, 28	5Fr, 1.5ml	10	R (GAMMA)
TWO-TAP HS CATHETER	5FR/7FR	28	5Fr, 1.5ml/7Fr, 3ml	10	R (GAMMA)



65 Parker St., Unit #2
 Newburyport, MA 01950
 Phone: (800) 531-2545

Fax: (978) 914-6106
 Email: info@kmmmedical.com

